




TinkerCell User Documentation

1. Biological components	2
Units	3
2. Connections	3
Visual details	4
3. Inserting components.....	5
4. Inserting connections	6
5. Secret Key Strokes.....	8
6. Analysis of models	9
7. View and edit model equations	9
8. View and edit parameters	11
Naming convention	11
Automatic correction	11
Viewing and changing parameters.....	12
9. Compartments	12
10. Events and forcing functions.....	13
11. Text based model construction	14
12. User home folder	14
13. Getting updates	15

14. Programming.....	15
Sharing code.....	15
Command-line	15
PySCeS	16

Biological components

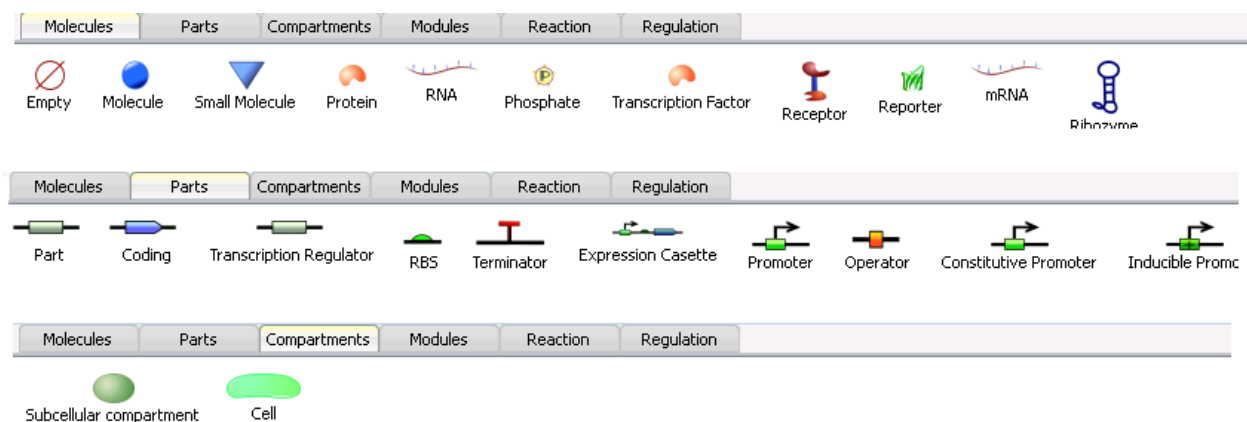

The catalog of biological components (shown above) are the fundamental objects that are used to construct models in TinkerCell. Each component has a unique set of attributes and a designated visual representation. The attributes characterize that part. For example, the attributes for a "Coding" part are its DNA sequence (a text string) and the length of the sequence (a number). The default visual representation of each part can be changed by the user.

The catalog is organized as a tree, although not displayed that way. The tree of components is organized into three major divisions: Molecules, Parts, and Compartments. The "empty" object is in a category by itself. The "Molecules" division represents components such as proteins, RNA, and metabolites. The "Parts" division represents components that belong on DNA, such as genes and promoters. The "Compartments" division represents items with volume that can contain other components. Compartments can include cells, sub-


cellular regions such as nucleus or mitochondrion, or other separated regions within the system.

The philosophy of TinkerCell is to make models well characterized. Well characterized models will enable TinkerCell to exchange information with biological databases in future. It also allows automation, as the "Modules" section later will demonstrate. To characterize models well, TinkerCell models also contains information about the type of item, such as "Promoter", "Protein", and "mRNA", instead of defining all components of a model as generic "variables" and "reactions". The list of components and their hierarchy are stored in an XML file called NodesTree.xml. Editing this file will alter the tree of components. This file will eventually be loaded from a database of known biological components. Users may click the right mouse button on components in the catalog to learn more about each item.

Units

Units are not defined for an entirely family of components (e.g. proteins) rather than individual items. The Settings menu at the top provides an option for setting the units for each type of component.

Connections

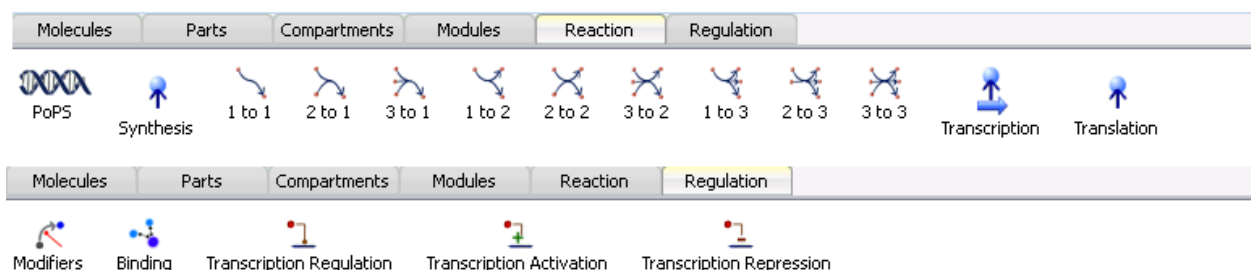


The catalog of connections (shown above) is the set of different types of connections that can be used to connect molecules and parts. Each type of connections can only be made between specific types of nodes. For example,

a reaction where one small molecule such as glucose is converted to another small molecule is classified as a "biochemical" reaction in TinkerCell. Sub-classes of "biochemical" reactions include the "1 to 1", "2 to 1", "2 to 2", and similar connections listed in the catalog. When a transcription factor regulates a promoter, it is classified under "transcription regulation", which can then be further classified as activation or repression. A transcription regulation cannot be used to connect two molecules together (TinkerCell will display an error). Some of these connections only allow two parts to be connected, while some allow different number of parts to be connected. The type and number of parts allowed for each connection can be viewed by clicking the right mouse button on one of the connections in the connections catalog.

Once again, the philosophy in TinkerCell is to make models well characterized. If all the reactions were simply called "reactions", it would be very difficult for TinkerCell to interact with databases or perform other tasks which requires knowledge about the type of connections.

Visual details

Connections can be decorated in many ways: the arrow heads can be changed using the "replace graphics" option in the toolbar; the "middle region" can be made visible and replaced using the "replace graphics" option again. The lines can be dashed, dotted, or solid. Modifier connections (bottom-most option in the connections tree) are used to connect a part to a connection. Below is an illustration.

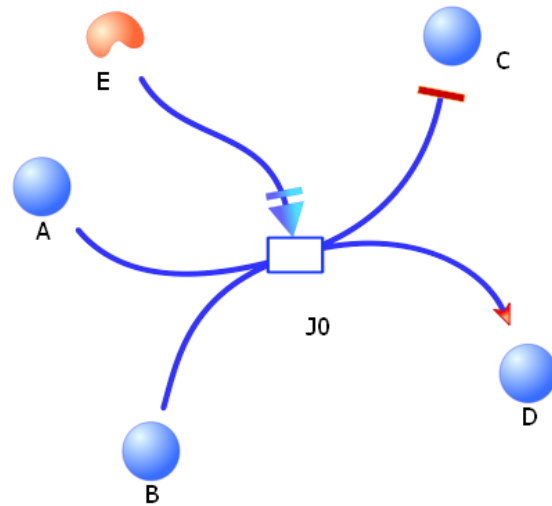

Inserting components

How to insert components:

Components are inserted by clicking on one of the buttons in the components catalog and then clicking on the canvas. You can also drag and drop items to the main canvas.

Extra details: what happens a component is inserted:

Various "plug-ins" in TinkerCell monitor the user's activity. When a component is inserted, the plug-ins carry out various operations. For example, one of the plug-ins called the "Parameters" retrieves information from the catalog of components and connections in order to define default *parameters* of each component. Some of the plug-ins automatically create reaction rate equations or insert sub-models.

Decorating components:

Using the "Add decorator" option in the context-menu (i.e. mouse right click or ctrl+mouse click on a Mac), it is possible to add visual annotations on components. These decorators will move along with the objects that they belong with. **Hold CTRL when moving objects to move them independently.**




--> using CTRL to move the decorators -->



Inserting connections

How to insert connections:

Connections are inserted in two ways. First, they can be inserted by dragging and dropping the connection on the canvas. Second, they can be inserted by first clicking on one of the connections and then selecting the components to connect, one by one. The connections tree imposes some restrictions about what type of components can be connected with the different types of connections. For example, a user may not connect a promoter and a gene using an "Allosteric inhibition" connection because the definition of "Allosteric inhibition" requires a small molecule and a protein.

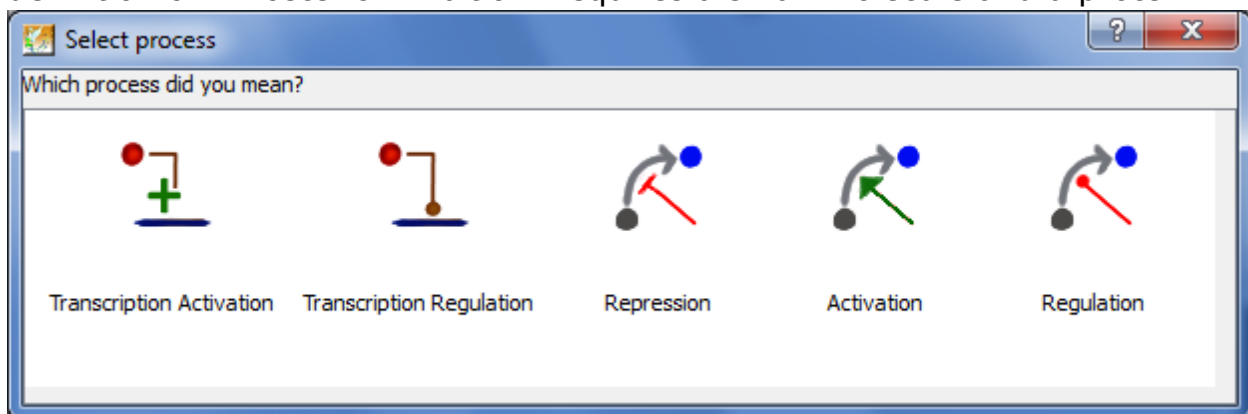


However, the user may connect any two components using generic connections such as "repression", "activation", or "1 to 1". For generic

connections, TinkerCell will ask the user whether they meant a more specific type of connection. For example, when the user connects a Transcription Factor to an Inducible Promoter, TinkerCell will ask which connection they mean, with the first option being the most specific option (in this case, Transcription Activation), as shown in the image to the right.

Process diagrams:

TinkerCell uses a method of modeling where a sub-models, or a process, can be represented by a single connection. When a user connects a small molecule to a protein using the "Allosteric inhibition" connection, then TinkerCell automatically loads a simple binding-unbinding type of model within that connection. The user can view or edit this sub-model by double-clicking (or pressing Enter) on the connection itself, as shown in the image below. Some of the components inside the sub-model will be the same as the components in the outside model. These components will be highlighted in gold color when the corresponding component is selected in the sub-model or the outside model.

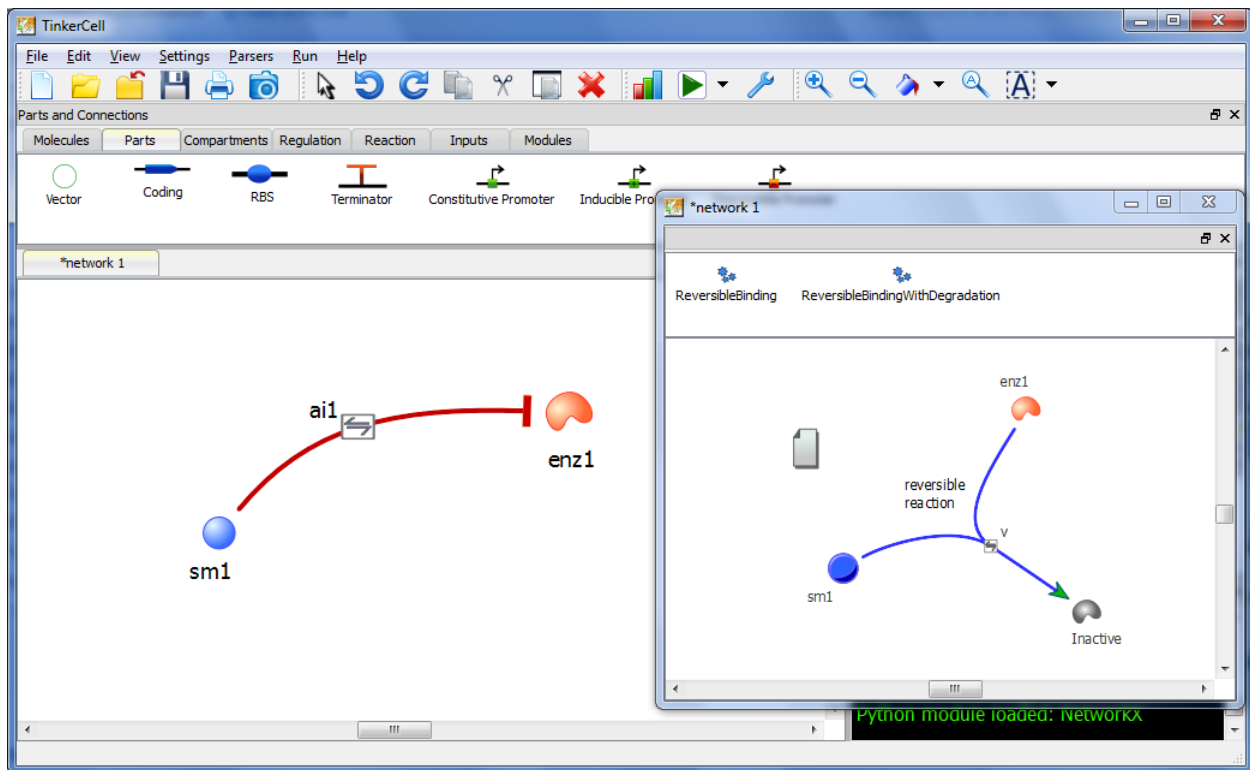

Secret Key Strokes










- Hold CTRL to move individual items such as a cell by itself without the contents within in.
- Hold CTRL and click on an empty part of the screen to drag the screen using the mouse
- Use PAGE-UP and PAGE-DOWN to change the ordering of selected items, i.e. send them back or bring them forward
- Use CTRL and double-click to add new control points on connections
- Hold CTRL and use the mouse wheel to zoom
- Use the + and - keys to zoom in and out
- Hold CTRL and the arrow keys to move the visible area

Analysis of models

TinkerCell provides numerous simulation and analysis functions. All the analysis functions are located in the "Programs" tab on the side (see image below). The same functions are also available in two other places: under the



icon on the toolbar and under the "Run" menu. Each function is listed under a category. For example, the "Simulate" category contains deterministic simulation and stochastic simulation. The "Steady state" category contains eigenvalue calculations, 2D and 1D steady state analysis, which can be used to study the input-output response of a system. Most of these simulation functions come with sliders, allowing users to see the affect each parameter has on the simulation result interactively.

-
- ▶  Simulation
 - ▶  Steady State
 - ▶  Metabolic Control Analysis
 - ▶  Network Structure
 - ▶  Optimization
 - ▶  Generate kinetics
 - ▶  Sequence
 - ▶  Database
 - ▶  Model Variations
-

View and edit model equations

There are several ways to examine and edit a model. The first method is to use the tools that are displayed when the user selects items on the screen

[illegible]

Initial Values	Parameters	Rates	
Initial values of selected items			
	name	value	fixed?
	s	1	floating
c.	r	1	floating
c.c2.	r	1	floating
c.c2.	p	1	floating
c.c2.	s	1	floating
c.	s	1	floating

Alternatively, the user may double-click or press enter on selected items to see the model summary window.

View and edit parameters

Parameters, also called "Numerical Attributes" in TinkerCell (only relevant for programmers), are constant values that may or may not be used in the model. The important fact is that parameters are local, meaning that they belong with an item. For example, reaction J0 might have parameter k0 and reaction J1 might also have a parameter called k0. The full name of each parameter is J0.k0 and J1.k0, respectively. When entering an equation, the full parameter name should be used. If the user enters "k0" by accident, TinkerCell will select one of the two: J0.k0 or J1.k0.

Naming convention


When an item has a parent item, such as a compartment, the parent's name is prepended to the name. For example, if part "A" is inside compartment or module "B", then the complete name of "A" is "B.A". If part A has a parameter named "k1", then the way to address this parameter in other parts of the model would be "B.A.k1".

Automatic correction

Every time the model is modified, TinkerCell generates an internal table with all the individual item names and parameter names as well as complete names. If a user enters an incomplete name, TinkerCell will search this table and use the closest matching complete name. If no such name exists, TinkerCell will create a new parameter with that name. The new parameter will belong with the item under consideration.

Viewing and changing parameters

The easiest way to view and edit any of the parameters is the "Model summary", shown to the right. This window lists all the items in the model and their respective parameters. Equations can also be edited, except for situations in which a reaction might have two or more rates (e.g. forward and backward), in which case this window will only show the first equation.

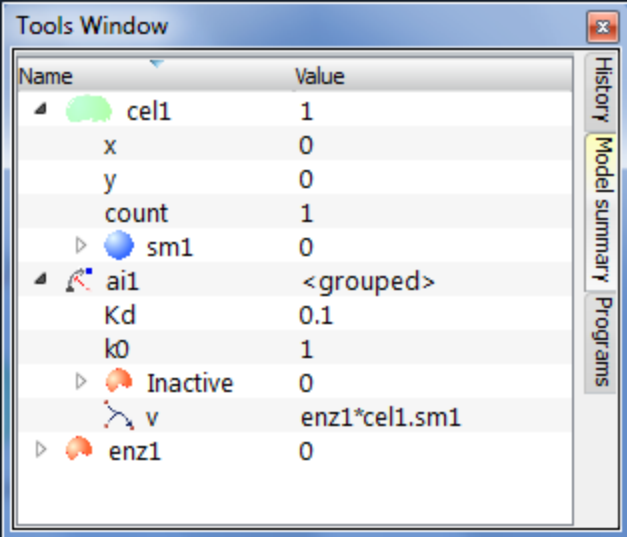


Name	Value
cel1	1
x	0
y	0
count	1
sm1	0
ai1	<grouped>
Kd	0.1
k0	1
Inactive	0
v	enz1*cel1.sm1
enz1	0

The main reason for using local parameters is because the parameters can be loaded from a database or from some other source which might store information about a biological part. If parameters did not "belong" to a part, it would be nearly impossible to determine which parameters of the system to modify when a single part is loaded from a database.

Compartments

Compartments represent a space (volume) inside of which certain items are localized. In the screenshot below, two compartments, c1 and c2 belong inside a larger compartment, the cell. While the names appear to be redundant, the "complete" name is unique, because the complete name has the parent item's name as a prefix, which is explained in the next section. The parent of an item is the compartment containing that item.

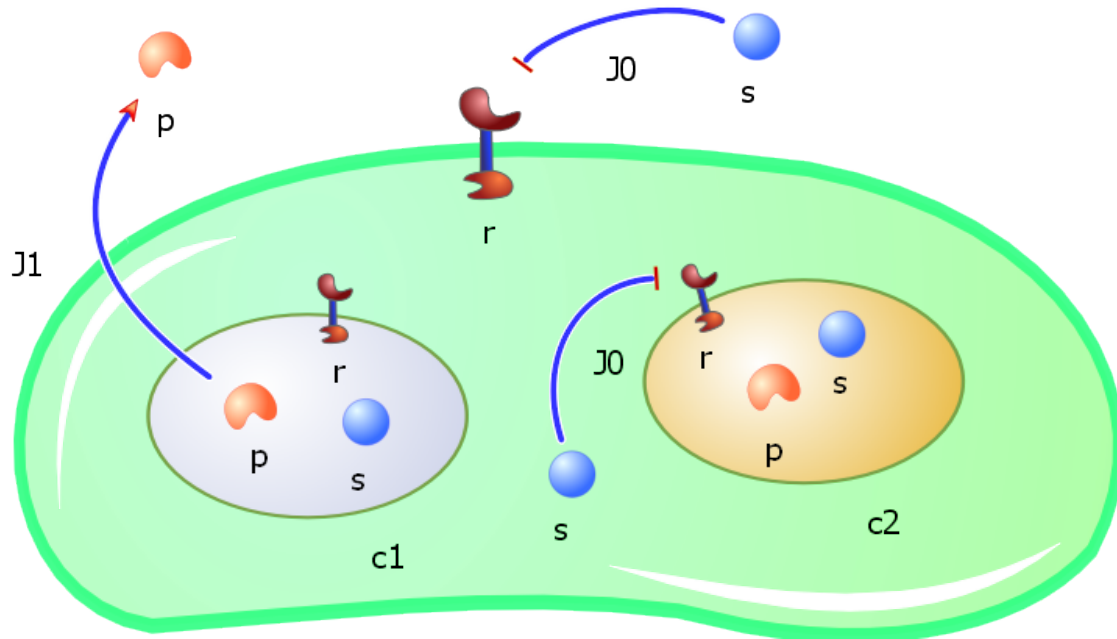

Events and forcing functions

There is an "Inputs" tab located next to the components and connections catalog. There are a list of forcing functions, such as sine wave, that are available here. There is also an icon for inserting events. Events are responses that are triggered under a specific condition. Insert an event on the screen and double-click in order to add or remove events.

Forcing functions can also be provided by selecting the "Functions" tool that is displayed when items are selected. This tool will list a set of functions defined inside the selected items. For example, if a user wants to assign a function to the molecule "s1", then the user should select this molecule, select the "Function" tool on the right, and then enter $s1 = s2 * (1 + \sin(\text{time}))$ or some other function of interest.

Text based model construction

TinkerCell supports graphical and text-based modeling. The text-based modeling is done via a third-party text parser (a plug-in). The text parser that is included with TinkerCell is called Antimony (antimony.sourceforge.net). Antimony is a modular model construction language. An example script is shown below.

Model M

```
E + S -> ES;   kf * E*S
```

```
ES -> E + S;   kb * E*S
```

```
ES -> EP;   k1 * ES
```

```
EP -> ES;   k2 * EP
```

```
EP -> E + P;   kcat * EP
```

```
kf = 1 + 0.6 *sin(time)
```

```
kb = 1.2
```

```
k1 = 0.1
```

```
k2 = 0.9
```

```
kcat = 0.3
```

```
end
```

User home folder

TinkerCell designated a folder as the user's home folder for temporary TinkerCell files. When TinkerCell loads plug-ins, it will also search this folder. So users without write-access to the main installation folder may use their TinkerCell home folder to place new plug-ins or Python scripts. By default, this folder is **<Documents folder>/TinkerCell**, where <Documents folder> is the default documents directory of the current user.

Getting updates

Use the Download button (down arrow) located in the toolbar to automatically download the latest set of modules and plugins that others have uploaded to the TinkerCell repository. The same option is also available in the Settings menu.

Programming

C, Python, and Octave functions can be written in TinkerCell, and these functions can be added to the list of programs in TinkerCell. See [TinkerCell C API](#). The programming language can be selected from the Settings menu. The programming window (wrench icon) can be used to write multi-line code. The console window can be used to write single-line (quick) code.

Sharing code

When code written in the coding window is saved, TinkerCell will provide an option to the user for uploading the code to a central repository. All other TinkerCell users will automatically get this piece of code when they start TinkerCell.

Command-line

The command-line, or console window, is used for three purposes: (1) it serves as a display window for errors and other information; (2) the user can enter Python or Octave commands in the console window; (3) the user can enter the name of any component or parameter in the model and see information about that component or value of that parameter.

PySCeS

All functions in the PySCeS module can be used inside TinkerCell's command prompt. PySCeS model is loaded by invoking **mod = tc2pysces.load()** (tc2pysces stands for TinkerCell-to-PySCeS). The returned variable **mod** will be a PySCeS model. Please see the [PySCeS Documentation](#) for all the functions available.